


# EGC Assessment: An introduction for Healthcare Providers

## Chest electrode and lead placements: Reference Sheet

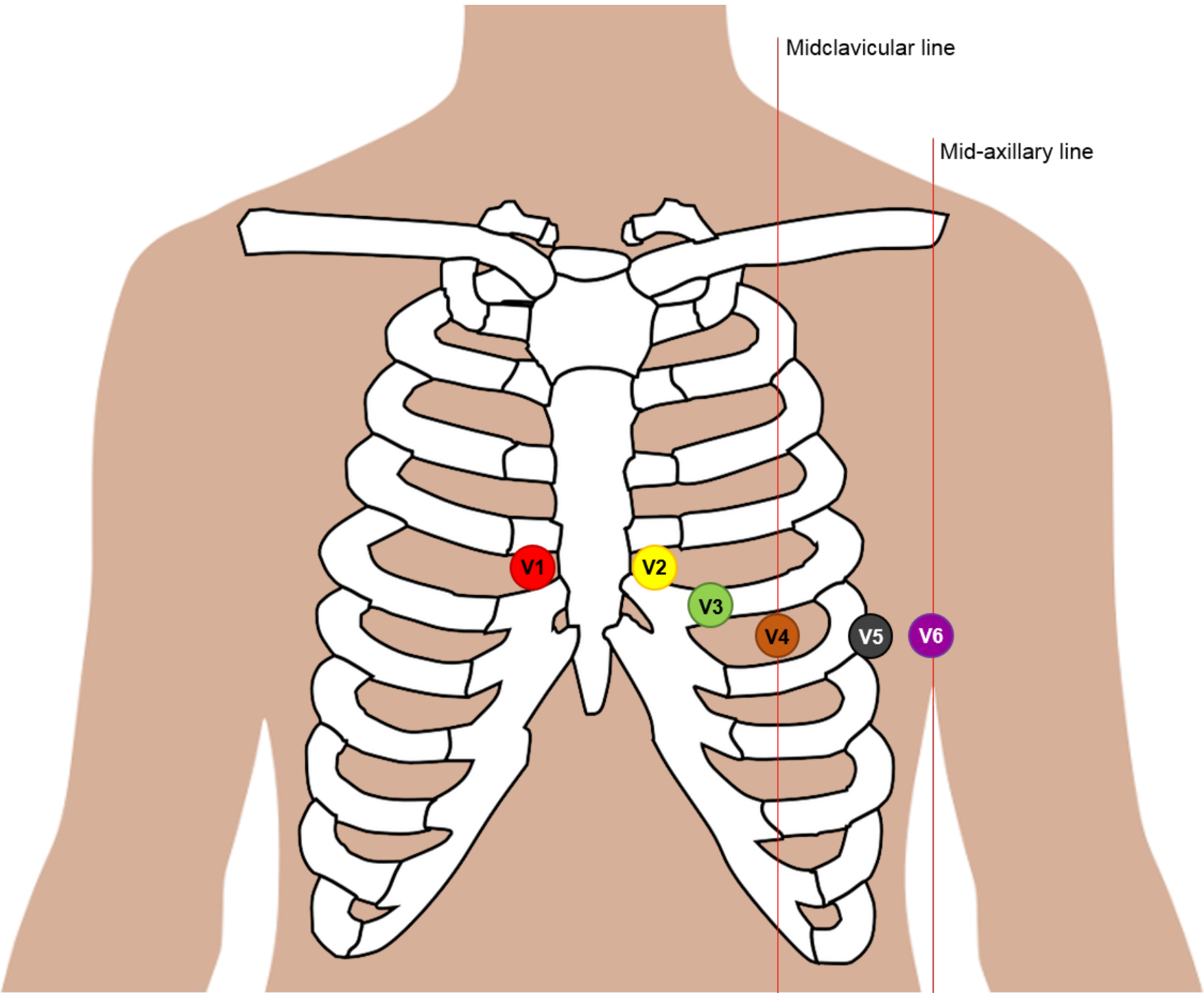

Label	Position of electrode / lead on body	Colour
V1	Fourth intercostal space at the right border of the sternum	Red
V2	Fourth intercostal space at the left border of the sternum	Yellow
V3	Midway between placement of V2 and V4	Green
V4	Fifth intercostal space at the midclavicular line	Brown
V5	Anterior axillary line on the same horizontal level as V4	Black
V6	Mid-axillary line on the same horizontal level as V4 and V5	Purple